


Sample Information

Name	:Tirzepatide		
Sequence	:Tyr-{Aib}-Glu-Gly-Thr-Phe-Thr-Ser-Asp-Tyr-Ser-Ile-{Aib}-Leu-Asp-Lys-Ile-Ala-Gln-{diacid-C20-gamma-Glu-(AEEA)2-Lys}-Ala-Phe-Val-Gln-Trp-Leu-Ile-Ala-Gly-Gly-Pro-Ser-Ser-Gly-Ala-Pro-Pro-Pro-Ser-NH2		
Pump A	:0.1%Trifluoroacetic in 100% water		
Pump B	:0.1%Trifluoroacetic in 100% acetonitrile		
Total Flow	:1ml/min		
Wavelength	:214nm		
Analytical column type	:Phenomenex Luna® 5 µm C18 (2) (4.6*250mm*5 µm)		
Dissolution method	:100%H2O		
Inj. Volume	: 30ul		
Time	Module	Action	Value
0.01	Pumps	B.Conc	45
25.00	Pumps	B.Conc	70
33.00	Pumps	B.Conc	100
38.00	Pumps	B.Conc	100
40.00	Pumps	B.Conc	45
50.00	Controller	Stop	

Chromatogram

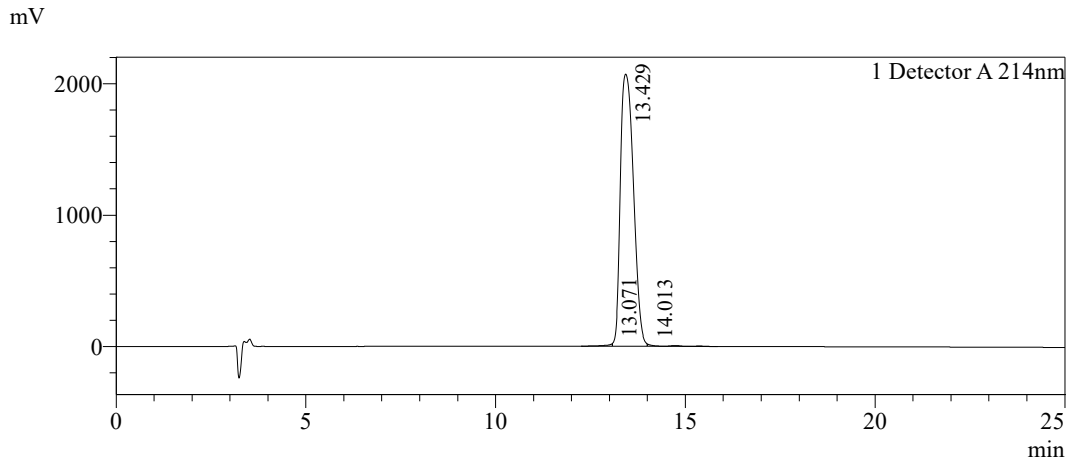

Peak Table

Detector A 214nm

Peak#	Ret. Time	Area	Height	Area%
1	13.071	223355	18979	0.446
2	13.429	49517862	2071039	98.886
3	14.013	334620	16853	0.668
Total		50075838	2106871	100.000