


Sample Information

Name :Acetyl Octapeptide-3  
 Sequence :Ac-EEMQRRAD-NH2  
 Pump A :0.1%Trifluoroacetic in 100% water  
 Pump B :0.1%Trifluoroacetic in 100% acetonitrile  
 Total Flow :1ml/min  
 Wavelength :214nm  
 Analytical column type :Phenomenex Luna® 5 µm C18 (2) (4.6\*250mm\*5 µm)  
 Dissolution method :100%H2O  
 Inj. Volume : 25ul

Time	Module	Action	Value
0.01	Pumps	B.Conc	5
25.00	Pumps	B.Conc	30
33.00	Pumps	B.Conc	100
38.00	Pumps	B.Conc	100
40.00	Pumps	B.Conc	5
50.00	Controller	Stop	

Chromatogram

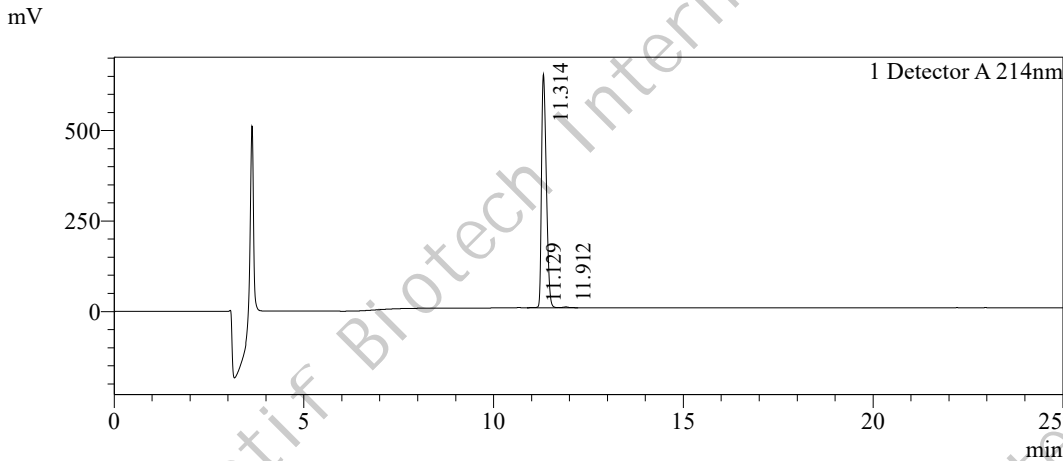

Peak Table

Detector A 214nm

Peak#	Ret. Time	Area	Height	Area%
1	11.129	11425	1245	0.199
2	11.314	5677395	645801	99.118
3	11.912	39078	2664	0.682
Total		5727898	649709	100.000

Email : marketing@motifbiotech.com

Web : https://www.motifbiotech.com

Address: B04, Building 17, No. 122, Yongan Rd, Suzhou National Hi-Tech District, PRC2