


Sample Information

Name :Palmitoyl Hexapeptide-12
 Sequence :Pal-VGVAPG
 Pump A :0.1%Trifluoroacetic in 100% water
 Pump B :0.1%Trifluoroacetic in 100% acetonitrile
 Total Flow :1ml/min
 Wavelength :214nm
 Analytical column type :Phenomenex Luna® 5 µm C18(2) (4.6*250mm*5 µm)
 Dissolution method :15%ACN+85%H2O
 Inj. Volume : 15ul

Time	Pumps	Action	Value
0.01	Module	B.Conc	75
25.00	Pumps	B.Conc	100
33.00	Pumps	B.Conc	100
38.00	Pumps	B.Conc	100
40.00	Pumps	B.Conc	75
50.00	Controller	Stop	

Chromatogram

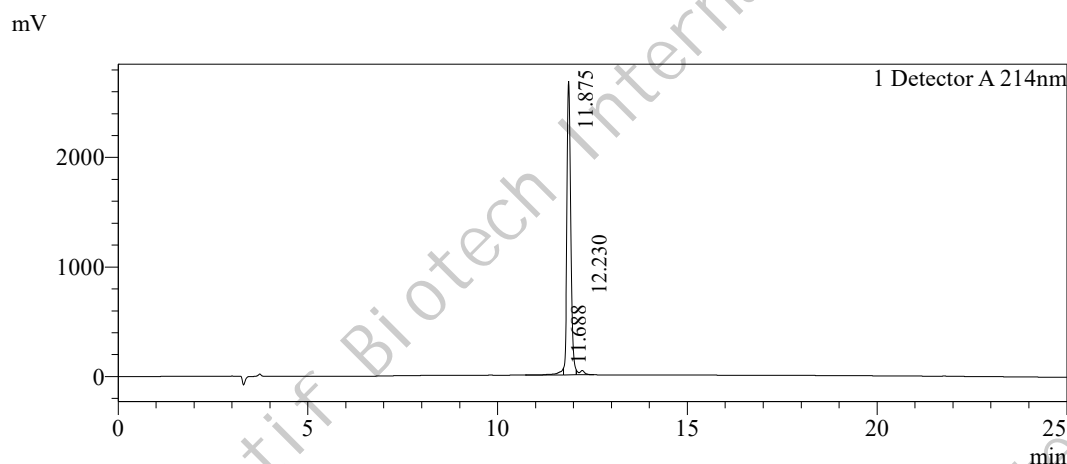

Peak Table

Detector A 214nm

Peak#	Ret. Time	Area	Height	Area%
1	11.688	479905	34180	2.390
2	11.875	19165650	2680725	95.443
3	12.230	435229	40244	2.167
Total		20080785	2755148	100.000

Web : <https://www.motifbiotech.com>

Email : marketing@motifbiotech.com

Add : B04, Building 17, No. 122, Yongan Rd, Suzhou National Hi-Tech District, PRC